

BREAST BIOPSY




Your partner in Biopsy and Special Needles

MEDONE ULTRA

PROGRAMMABLE AUTOMATIC DISPOSABLE BIOPSY SYSTEM

- 25 mm penetration depth, 20 mm sample notch
- Automatic or sequential firing modes
- Echogenic tip for ultrasound guided procedures
- Powerful spring firing mechanism
- Centimeter markings

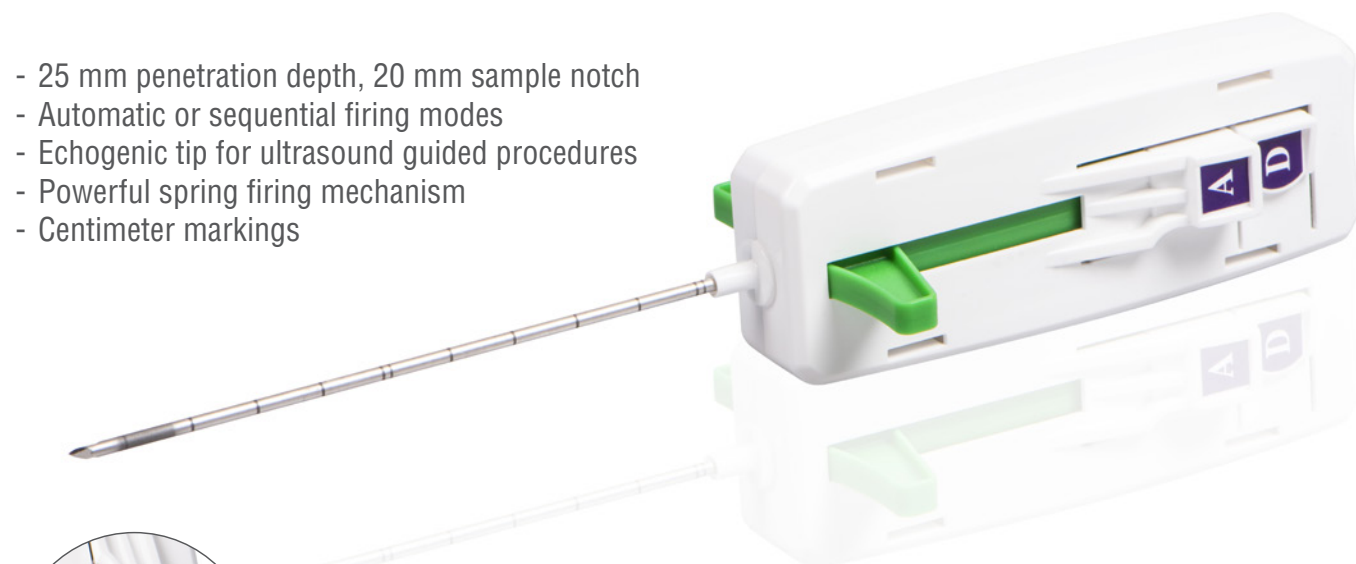

Dual firing mode:

- D** Delay firing mode to visualize needle stylet sample notch placement within a lesion before activating the cutting cannula
- A** Automatic firing mode release the stylet and cannula in rapid sequence capturing specimen by pushing a single button

Coaxial introducer for performing multiple biopsy procedures including Blunt Tip Stylet can reduce the risk of damage to vasculature or other vital structures



GAUGE	COLOR	LENGTH (MM)		
12	Light blue	12 100	12 160	12 200
14	Green	14 100	14 160	14 200
16	White	16 100	16 160	16 200
18	Pink	18 100	18 160	18 200
20	Yellow	20 100	20 160	20 200

Other lengths available as special orders;
minimum quantities may apply


MEDONE

AUTOMATIC DISPOSABLE BIOPSY SYSTEM

- One hand usage
- 22mm penetration depth, 19mm sample notch
- Dual firing triggers in proximal and distal positions
- Echogenic tip for ultrasound guided procedures
- Centimeter markings
- Can be used with coaxial introducer needle
- One-hand use leaves the physician able to operate the ultrasound probe
- Ergonomic shape and optimal position firing buttons create a user-friendly device



ONE HAND OPERATION

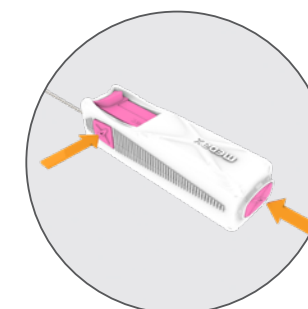

OPTIONAL FIRING TRIGGERS
POSITIONED ON SIDE AND REAR


GAUGE	COLOR CODE	LENGTH (MM)				
12	Light blue	12 100	12 160			
14	Green	14 100	14 160	14 200		
16	White	16 100	16 160	16 200		
18	Pink	18 100	18 160	18 200	18 250	18 300
20	Yellow	20 100	20 160	20 200		

Other lengths available as special orders; minimum quantities may apply


MEDCUT

AUTOMATIC DISPOSABLE BIOPSY SYSTEM

- 22mm penetration depth, 19mm sample notch
- Firing trigger in proximal position
- Extra sharp tip for safe insertion and penetration
- Echogenic tip for ultrasound guided procedures
- Centimeter markings

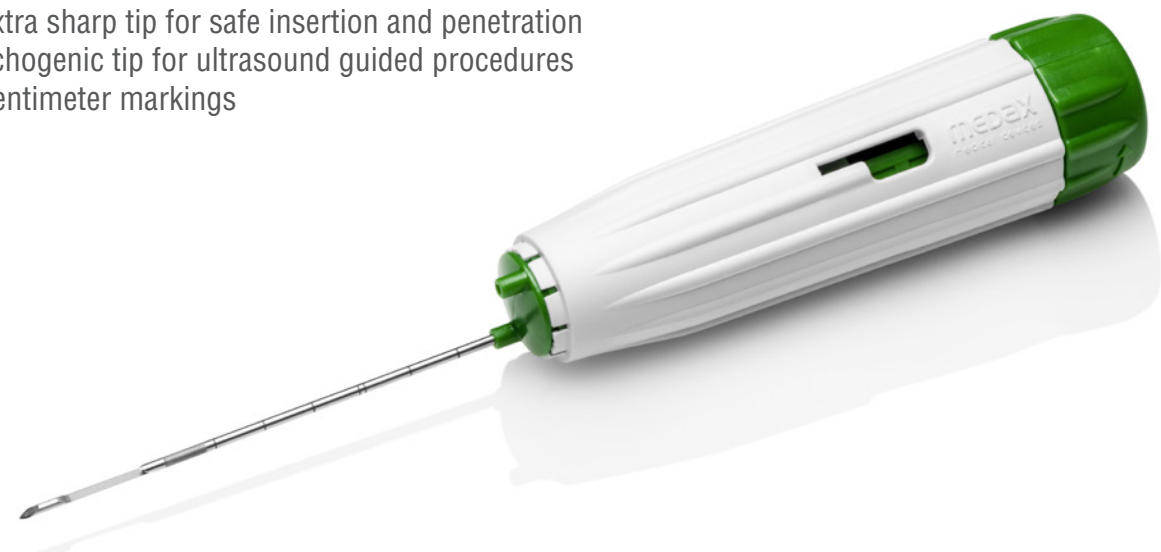

Side needle position in order to make Medcut compatible with any kind of abdominal and transrectal probe

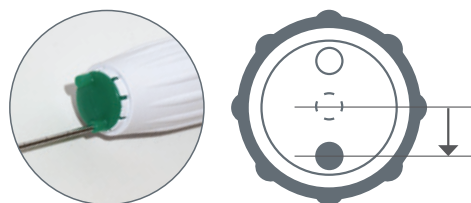


GAUGE	COLOR CODE	LENGTH (MM)				
14	● Green	14 100	14 160			
16	○ White	16 100	16 160	16 200		
18	● Pink	18 100	18 160	18 200	18 250	18 300
20	● Yellow	20 100	20 160	20 200		

Other lengths and tips available as special orders; minimum quantities may apply

MEDEASY

SEMI-AUTOMATIC DISPOSABLE BIOPSY SYSTEM

- Ergonomic lightweight design with open finger rings
- Adjustable penetration depth (10mm or 20mm)
- Echogenic treatment for ultrasound guided procedures
- Ideal for CT guided procedures.
- Centimeter markings
- Unique handle design provides multiple ways to grip the device
- 10mm or 20mm penetration depth indicator displays on the handle

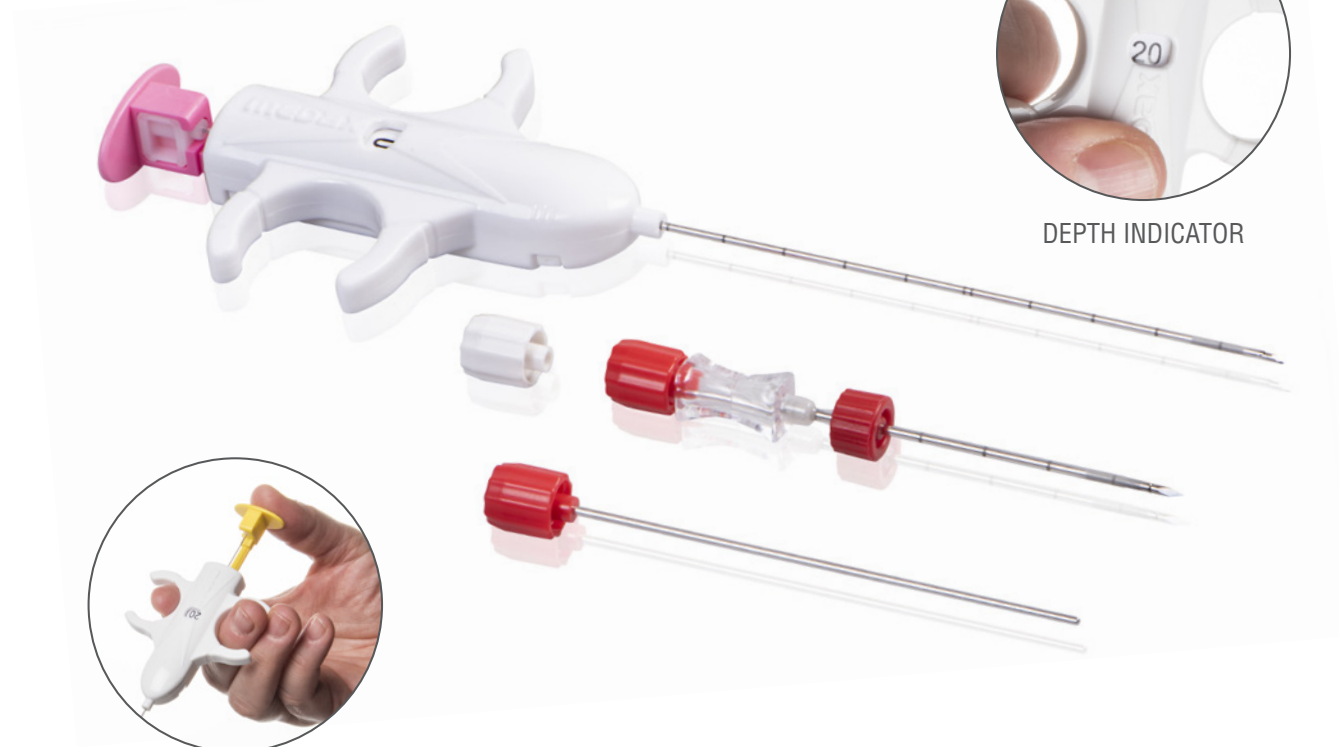

DEPTH INDICATOR

ERGONOMIC HANDLE

GAUGE	COLOR CODE	LENGTH (MM)				
12	● Light blue		12 100	12 160	12 200	
14	● Green		14 100	14 160	14 200	
16	○ White	16 080	16 100	16 160	16 200	
18	● Pink	18 080	18 100	18 160	18 200	18 250
20	● Yellow	20 080	20 100	20 160	20 200	20 250

Other lengths and tips available as special orders; minimum quantities may apply.

MEDWIRE

REPOSITIONABLE WIRE LOCALIZATION NEEDLE FOR NON-PALPABLE BREAST LESIONS

- Easy wire release
- Extremely sharp tip for atraumatic penetration
- Echogenic treatment for ultrasound guided procedures
- Centimeter markings

O-TYPE HOOK WIRE

J-TYPE HOOK WIRE

	GAUGE	Ø WIRE (MM)	LENGTH (MM)			
O-Type Hook Wire	20	0,40	20 060	20 090	20 120	
J-Type Hook Wire	20	0,40	20 030	20 050	20 070	20 100

Other lengths available as special orders; minimum order quantities may apply

NON-REPOSITIONABLE WIRE LOCALIZATION NEEDLE FOR NON-PALPABLE BREAST LESIONS

- Extremely sharp tip for atraumatic penetration
- Echogenic treatment for ultrasound guided procedures
- Centimeter markings and sliding stopper
- 5 year shelf life

SINGLE HOOK

STIFFENED VERSION

MARKED VERSION

VERSION HOOK	GAUGE Ø	WIRE (MM)	LENGTH (MM)		
L-Single	20	0,25-0,30	20 080	20 100	20 120

Other lengths and wire configuration available as special orders; minimum quantities may apply

ORDERING INFORMATION:

MEDONE ULTRA

M	U	GAUGE	LENGTH		
				0	0

INTRODUCER COAXIAL NEEDLE FOR MEDONE ULTRA
WITH BLUNT TIP STYLET

C	T	GAUGE	LENGTH		
				1	0

MEDONE ULTRA WITH COAXIAL NEEDLE INCLUDING
BLUNT TIP STYLET

M	U	GAUGE	LENGTH		
				K	0

MEDONE

M	O	GAUGE	LENGTH		
				0	0

INTRODUCER COAXIAL NEEDLE FOR MEDONE
WITH BLUNT TIP STYLET

C	O	GAUGE	LENGTH		
				1	0

MEDONE WITH COAXIAL NEEDLE
INCLUDING BLUNT TIP STYLET

M	O	GAUGE	LENGTH		
				K	0

MEDCUT

M	C	GAUGE	LENGTH		
				0	0

INTRODUCER COAXIAL NEEDLE FOR MEDCUT
WITH BLUNT TIP STYLET

C	M	GAUGE	LENGTH		
				1	0

MEDCUT WITH COAXIAL NEEDLE
INCLUDING BLUNT TIP STYLET

M	C	GAUGE	LENGTH		
				K	0

MEDEASY WITH COAXIAL NEEDLE
INCLUDING BLUNT TIP STYLET

M	L	GAUGE	LENGTH		
				K	0

MEDWIRE O-TYPE HOOK WIRE

X	W	GAUGE	LENGTH		
				0	0

MEDWIRE J-TYPE HOOK WIRE

X	W	GAUGE	LENGTH		
				0	J

MEDWIRE SINGLE HOOK WIRE

X	W	GAUGE	LENGTH		
				0	S

MEDWIRE SINGLE HOOK WIRE MARKED VERSION

X	W	GAUGE	LENGTH		
				1	S

MEDWIRE SINGLE HOOK WIRE STIFFENED VERSION

X	W	GAUGE	LENGTH		
				2	S



MEDAX S.R.L. Unipersonale

Via Sandro Pertini, 4 • 41039 • San Possidonio (MO) • Italy

Company direct No. : +39 0535 1812757 • customercare@medax.it • www.medax.it • www.medax-usa.com

Vat N. /Fiscal Code N. Iscriz. Reg. Impr. MN 02669860369 • N. REA: MO 403036 • Capitale Sociale Euro 100.011,00 i.v.



In accordance with the requirements of the medical device 93/42/EEC directive and its relevant updates.
All products undergo intensive clinical testing and are fully EC and FDA approved.