

CT/US BIOPSY



Your partner in Biopsy and Special Needles

MEDEASY
SEMI-AUTOMATIC DISPOSABLE BIOPSY SYSTEM

Semi-automatic spring loaded biopsy system designed to cover a wide range of soft tissue core biopsy procedures.

Key feautures:

- Ergonomic lightweight design with open finger rings
- Adjustable penetration depth (10mm or 20mm)
- Echogenic treatment for ultrasound guided procedures
- Centimeter markings
- Can be used with coaxial introducer needle (included)
- 5 year shelf life

Additional feautures:

- Ideal for CT guided procedures.
- Unique handle design provides multiple ways to grip the device
- 10mm or 20mm penetration depth indicator displays on the handle
- Blunt Tip Stylet included with coaxial introducer can reduce the risk of damage to vasculature or other vital structures
- Adapter for coaxial cannula assures accurate biopsy needle positioning when using the 10mm penetration depth



COAXIAL BLUNT TIP STYLET



ERGONOMIC HANDLE



DEPTH INDICATOR

GAUGE	COLOR CODE	LENGTH (MM)				
14	● Green		14 100	14 160	14 200	
16	○ White	16 080	16 100	16 160	16 200	
18	● Pink	18 080	18 100	18 160	18 200	18 250
20	● Yellow	20 080	20 100	20 160	20 200	20 250

Other lengths and tips available as special orders; minimum quantities may apply.

Ordering information:

MEDEASY WITH COAXIAL NEEDLE INCLUDING BLUNT TIP STYLET

Single box: 10 PCS.

M L GAUGE LENGTH K O

MEDONE
AUTOMATIC DISPOSABLE BIOPSY SYSTEM

Automatic, disposable biopsy system designed to cover a wide range of soft tissue core biopsy procedures

Key features:

- One hand usage
- 22mm penetration depth, 19mm sample notch
- Dual firing triggers in proximal and distal positions
- Echogenic tip for ultrasound guided procedures
- Centimeter markings
- Can be used with coaxial introducer needle (included)
- 5 year shelf life



GAUGE	COLOR CODE	LENGTH (MM)			
14	● Green	14 100	14 160		
16	○ White	16 100	16 160	16 200	
18	● Pink	18 100	18 160	18 200	18 250
20	● Yellow	20 100	20 160	20 200	

Other lengths available as special orders; minimum quantities may apply

Additional features:


- Powerful spring firing mechanism obtains optimal quality biopsy samples
- Extra sharp stylet tip and cutting edge cannula for precise penetration and high quality sample collection
- Ergonomic shape and optimal position firing buttons create a user-friendly device
- One-hand use leaves the physician able to operate the ultrasound probe
- Lightweight, compact design, when used with optional coaxial introducer needle, makes the system ideal for CT guided procedures as well as performing multiple biopsy procedures on the same site with minimal trauma
- Blunt tip stylet with coaxial introducer can reduce the risk of damage to vasculature or other vital structures




One Hand operation



Intuitive, user-friendly
arming mechanism

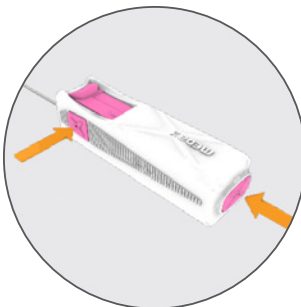

Optional firing triggers
positioned on side and rear

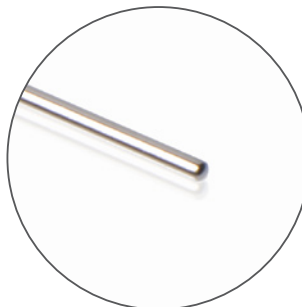

Coaxial Blunt Tip Stylet

Ordering information:

MEDONE

GAUGE		LENGTH			
M	O			0	0

MEDONE WITH COAXIAL NEEDLE INCLUDING BLUNT TIP STYLET

GAUGE		LENGTH			
M	O			K	0


Single Box: 10 PCS

MEDONE ULTRA

PROGRAMMABLE AUTOMATIC DISPOSABLE BIOPSY SYSTEM

Key feautres:

- One hand usage
- 25 mm penetration depth, 20 mm sample notch
- Automatic or sequential firing modes
- Echogenic treatment for ultrasound guided procedures
- Centimeter markings
- Can be used with coaxial introducer needle (included)
- 5 year shelf life

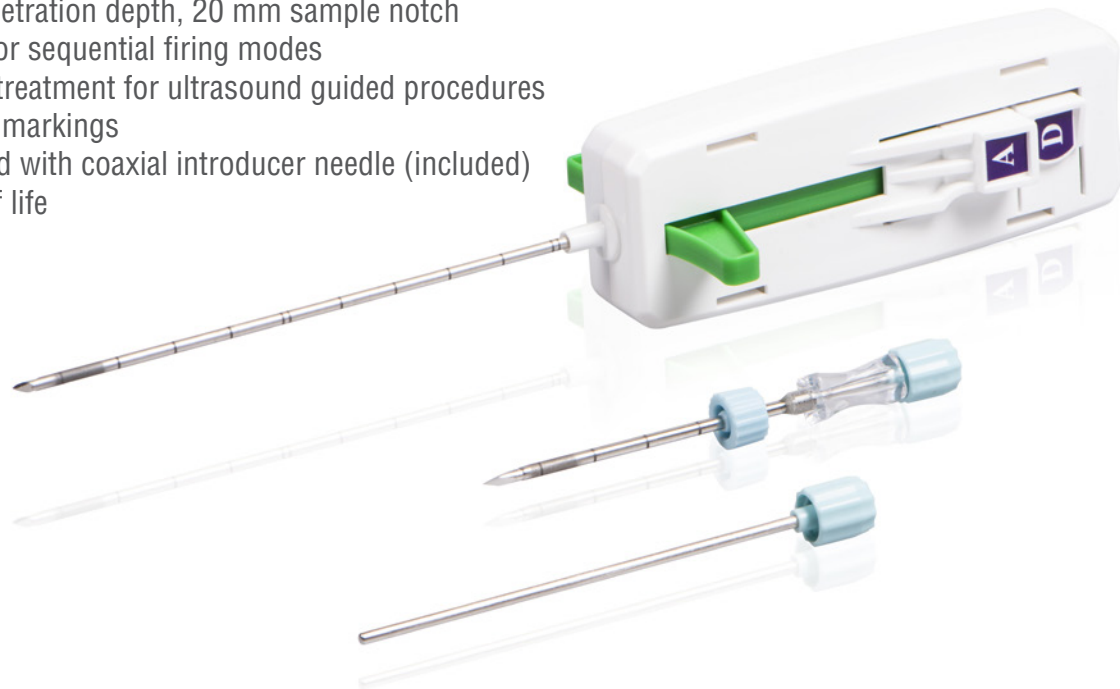

GAUGE	COLOR	LENGTH (MM)				
12	Light blue	12 100	12 120	12 160		
14	Green	14 100	14 120	14 160		
16	White	16 100		16 160	16 200	
18	Pink	18 100		18 160	18 200	18 250
20	Yellow	20 100		20 160	20 200	

Other lengths available as special orders;
minimum quantities may apply

Ordering information:

MEDONE ULTRA WITH COAXIAL NEEDLE
INCLUDING BLUNT TIP STYLET

MU

GAUGE

LENGTH

K0

Single Box: 5 PCS

FNA

FINE NEEDLE ASPIRATION FOR CYTOLOGICAL
AND MICROHISTOLOGICAL BIOPSY

Fine needle aspiration for cytological and
microhistological biopsy with different tip versions:
CHIBA, WESTCOTT and FRANSEEN

Key feautres:

- Transparent Luer Lock hub
- Sharpened cannula for safe and painless insertion
and penetration
- Internal and external echogenic treatment for correct
positioning under ultrasound guidance
- Centimeter markings and sliding depth stop
- 5 year shelf life

CHIBA needle:

Disposable needle used for
cytological biopsy or local
anesthesia administration

WESTCOTT and FRANSEEN
needles:

Disposable needle used for
cytological/microhistological
biopsies

GAUGE	COLOUR	LENGTH (MM)					
18	Pink	18 100	18 160		18 200		
20	Yellow	20 100	20 160		20 200		
22	Black	22 100	22 160	22 180	22 200	22 250	

Other gauge sizes and lengths available as special
orders; minimum quantities may apply



Ordering information:

CHIBA NEEDLE WITH ECHOGENIC TREATMENT

C

H

GAUGE

LENGTH

01

WESTCOTT NEEDLE WITH ECHOGENIC TREATMENT

W

E

GAUGE

LENGTH

01

FRANSEEN NEEDLE WITH ECHOGENIC TREATMENT

F

R

GAUGE

LENGTH

01

Single box: 20 PCS.

**MEDAX S.R.L. Unipersonale**

Via Sandro Pertini, 4 • 41039 • San Possidonio (MO) • Italy

Company direct No. : +39 0535 1812757 • customercare@medax.it • www.medax.it • www.medax-usa.com

Vat N. /Fiscal Code N. Iscriz. Reg. Impr. MN 02669860369 • N. REA: MO 403036 • Capitale Sociale Euro 100.011,00 i.v.



In accordance with the requirements of the medical device 93/42/EEC directive and its relevant updates.
All products undergo intensive clinical testing and are fully EC and FDA approved.